

## ARCADIS Avantic

**SP**

### Installation Instructions

System

Unpacking the Basic Unit

© Siemens AG 2005

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

**Document revision level**

The document corresponds to the version/revision level effective at the time of system delivery. Revisions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

**Disclaimer**

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such services.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.

1

Unpacking the Basic Unit

4

Work aids . . . . . 4

Workflow . . . . . 5

Changes to Previous Version . . . . . 11

## Work aids

Required tools:

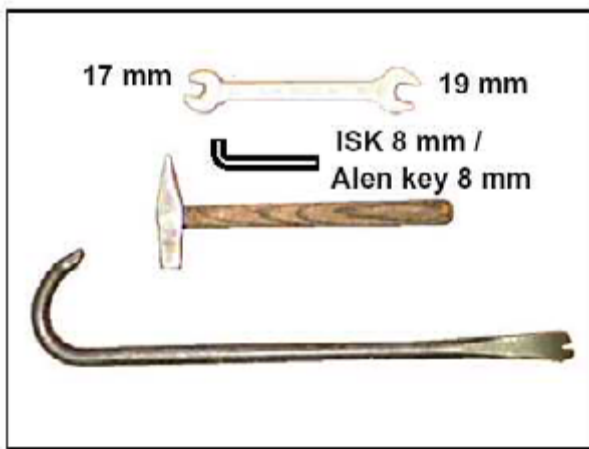


Fig. 1: Tools\_

## Workflow

- Cut the fastening strap, remove the ramp from the palett (1/Fig. 2 / p. 5), and lift out the board (2/Fig. 2 / p. 5).
- Unscrew the support (3/Fig. 2 / p. 5) and remove it from the palett.

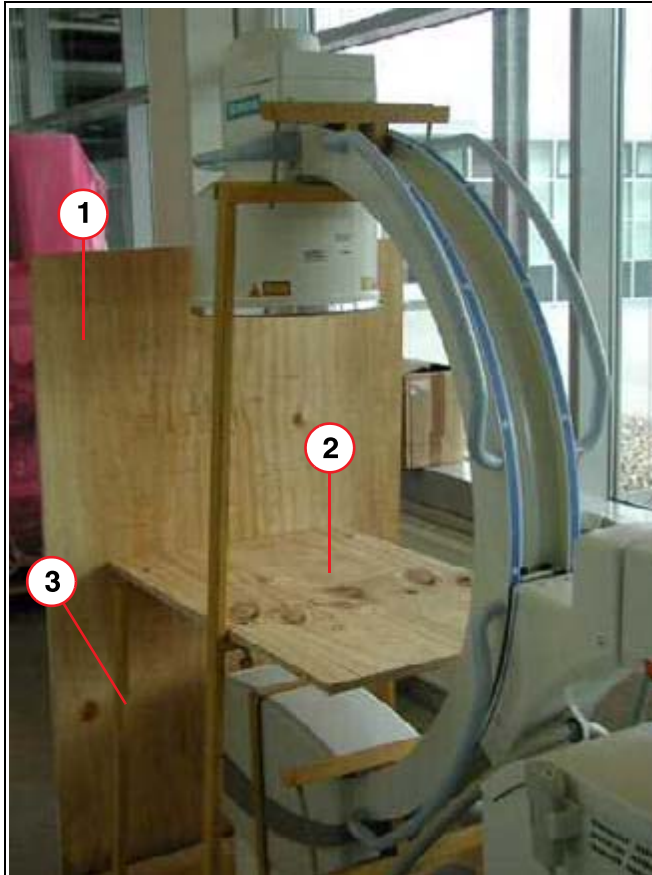


Fig. 2: Securing\_

- Unscrew and remove the retaining clamp (1/Fig. 3 / p. 6).



Fig. 3: \_Retaining clamp

- Unscrew the holding block (1/Fig. 4 / p. 6).



Fig. 4: \_Holding block

- Unscrew the retaining clamp (2/Fig. 5 / p. 7) on the left and right (screws).

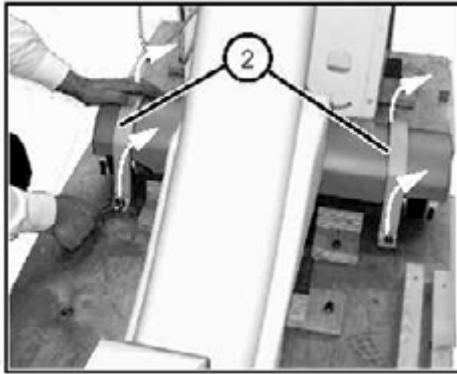


Fig. 5: Retaining clamp\_2

- Unscrew and remove the transport clamp (1/Fig. 6 / p. 7).

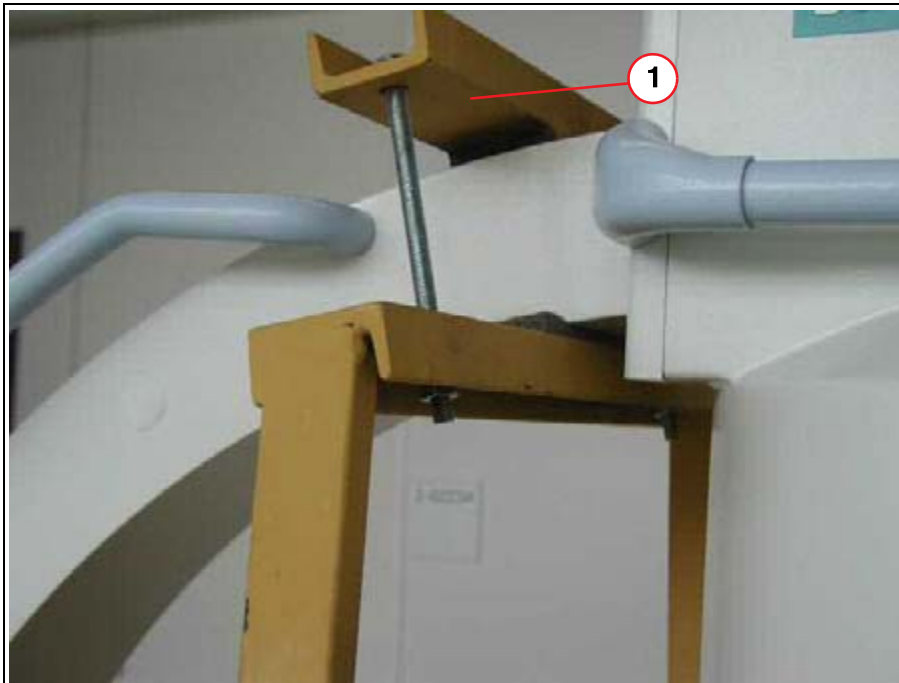


Fig. 6: \_Transport clamp

- Unscrew and remove the retaining clamp (1/Fig. 7 / p. 8) from the C-arm. Unscrew the support (2/Fig. 7 / p. 8).



Fig. 7: Support



- Unscrew and remove the retaining clamp (1/Fig. 8 / p. 9) (Styrofoam, (2/Fig. 8 / p. 9)).

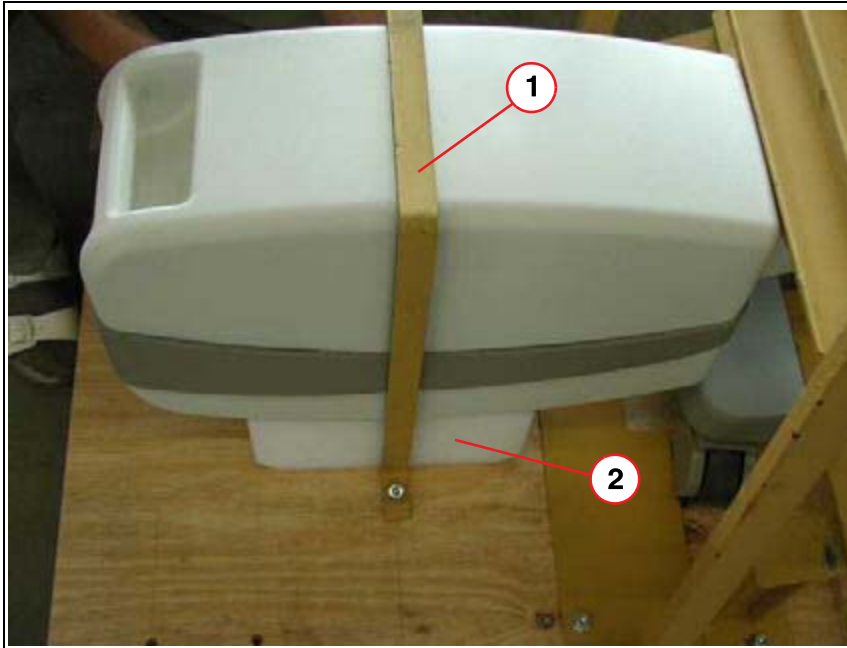


Fig. 8: Retaining clamp\_tube assembly

- Carefully pull the support toward the front (Fig. 9 / p. 9).



Fig. 9: Support

- Unscrew and remove the transport clamp (1/Fig. 10 / p. 10).

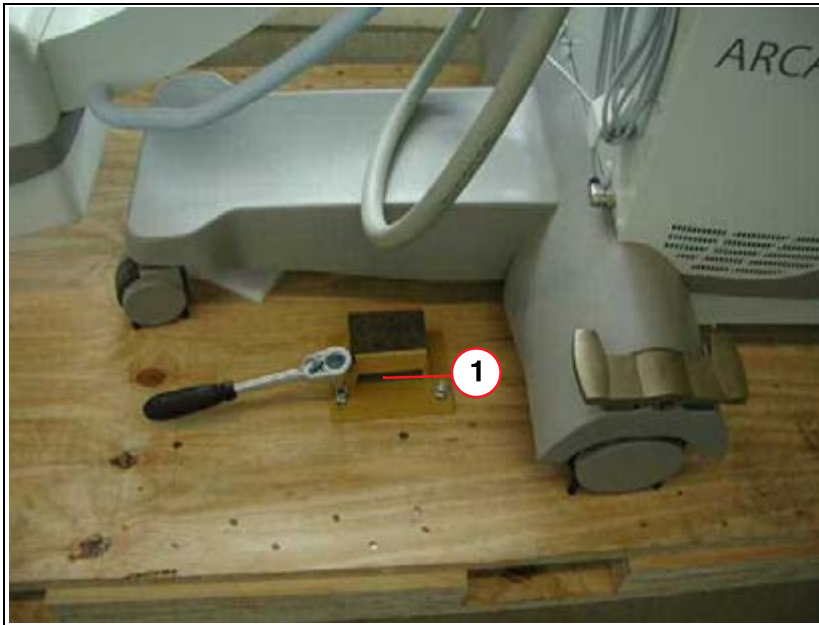


Fig. 10: Transport clamp\_

- Nail the ramp (1/Fig. 11 / p. 10) to the pallet.
- Release the foot brake and pull the unit in the direction of the pallet ramp .



Fig. 11: Ramp

## Changes to Previous Version

n.a.; initial version

